

## Electrothermal Slide Dry Bench




This Slide Drying Bench aids the preparation of microscope slides at the specimen mounting stage. Accepting up to 50 slides, it has the facility for drying slides in different orientations.

### Capacity

- Up to 50 slides
- The slides can be arranged on the bench in a number of ways as illustrated.

- Across drying bars
- Angled from bars
- Flat without bars

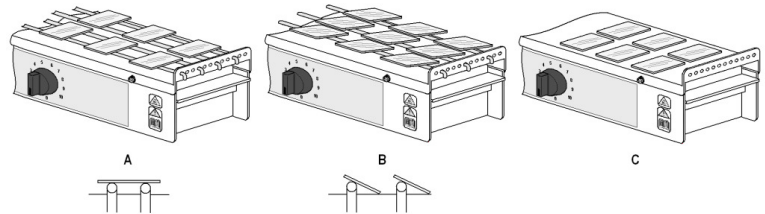

### Temperature

- Up to 100°C

### Pilot Light

- Clear to indicate “Power On” and Amber to indicate “Heater On”

Built-in energy regulator controlling temperature up to 100C.

Replacement elements easy fitted by user

**Warranty:** One year, parts and labour

**Certifications:** CE

Specifications		
<b>Material</b>	Aluminium Base and Top Stove Painted	
<b>Max. Element Temperature</b>	Up to 100°C	
<b>Cat. No.</b>	<b>Capacity</b>	<b>Electrical Requirement</b>
MH6616	50 slides	230V 50/60Hz
MH6616X1	50 slides	115v 50/60Hz
MH6616X6	50 slides	230v 50/60Hz, EU Plug